

syngo MammoReport

SP

CAD-Options

© Siemens AG 2005

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

#

Print No.: SPB7-250.812.34.01.02#

Replaces: n. a.

English

Doc. Gen. Date: 09.05

Document revision level

The document corresponds to the version/revision level effective at the time of system delivery. Revisions of hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

Purpose of this document

This manual is intended to provide technical service departments in subsidiaries and regional units with the information required to enable smooth start-up and maintenance of CAD-products.

Using CAD-option with Novation and MammoReport

CAD-license at MammoReport

Using the CAD-option requires a CAD-license at the MammoReport, depending on the software-version of MammoReport. This is for version VA20: SW-option CAD (8646023) supporting only R2-CAD, for version VB10 and higher: SW-option CAD V3.0 (10093545) supporting R2-CAD and iCAD.

CAD-workflow

Up to four Novations can deliver data for processing to the CAD-option, which is in case of R2 a UNIX-based workstation and in case of iCAD a Windows2000-based workstation.

It takes about one minute to process one case. The CAD-workstation sends the CAD-overlay to MammoReport and, if available, to the PACS-store.

Protection

R2-CAD functionality is protected by a USB-dongle that looks similar to a USB-memory-stick.

iCAD functionality is protected by license-strings, one license string is always available at the workstation to work as minimum configuration with one Novation, further (up to four) Novations require additional license-strings.

IP-addresses

Both CAD-workstations are delivered with a default IP-address and Subnet Mask that has to be modified to the IP-address and Subnet Mask given by the network administrator. Connection to the CAD-workstations is established via a net-interface, or if other servers share a KVP-switch by usage of that KVP-switch, or by adding a keyboard, mouse and SVGA-monitor to the R2-CAD or adding a SVGA-monitor to the iCAD (iCAD is delivered with keyboard and mouse).


Configuration

To configure the CAD-option, the AE-Title (case-sensitive!), IP-addresses and Subnet Mask of the Novation(s), the MammoReport(s) and, if available, the PACS-store(s) are required.

Additional the device-ID of the Novation(s) (serial number of Novation system) is needed to configure the CAD-option. The serial number of the Novation system can be obtained from the input in the syngo configuration page "site information":serial number".

In sequence the Novation(s), the MammoReport(s) and, if available, the PACS-store(s) need to be configured with the corresponding AE-Titles (case-sensitive!) IP-addresses and Subnet Masks according to the planned clinical workflow.

Port usage

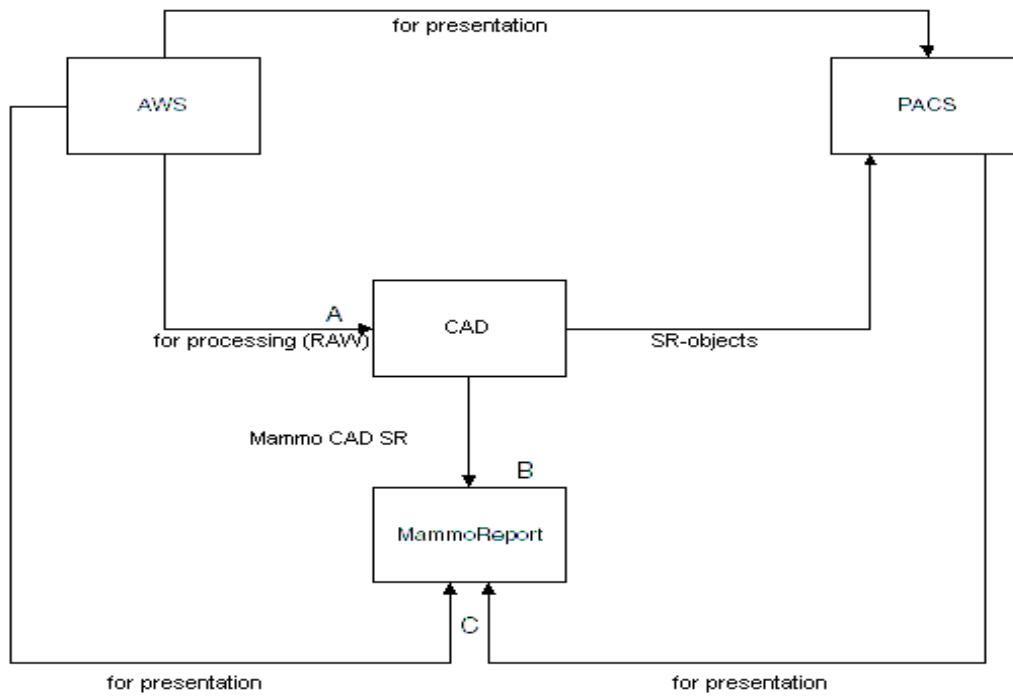

		R2	iCAD
	Objects	AET / Port	AET / Port
A	Image (Raw)	DICOMINBOX / 1115	AE_TITLE_CAD / 104
B (SW VA20)	Image + SR	5104	n. a.
B (SW \geq VB10)	SR (only!)	SCR_SYNGO / 50104	SCR_SYNGO / 50104
C	Image	5104	5104

Tab. 1

Changes since the previous version

Initial version



This page intenionally left blank